

How to install and setup the Bluetooth connection of g.MOBllab+ for paired Bluetooth devices, Version 3.10.01


The purpose of this tutorial is to instruct the user to set up and configure the paired Bluetooth connection to g.MOBllab+ based on the Ezurio Bluetooth USB dongle. This applies to Microsoft Windows 7 Professional English Bluetooth driver stacks. The device pairing has only to be done when the device is the first time installed or anytime you uninstalled g.MOBllab+.

1. Insert the Bluetooth dongle into a free USB slot. Wait until the Microsoft Bluetooth driver is installed automatically and the Bluetooth icon appears in the Notification area in the bottom right corner of your screen.



Double-click this icon.


2. In the **Bluetooth Devices** dialog click on the **Add...** button on the **Devices** tab.

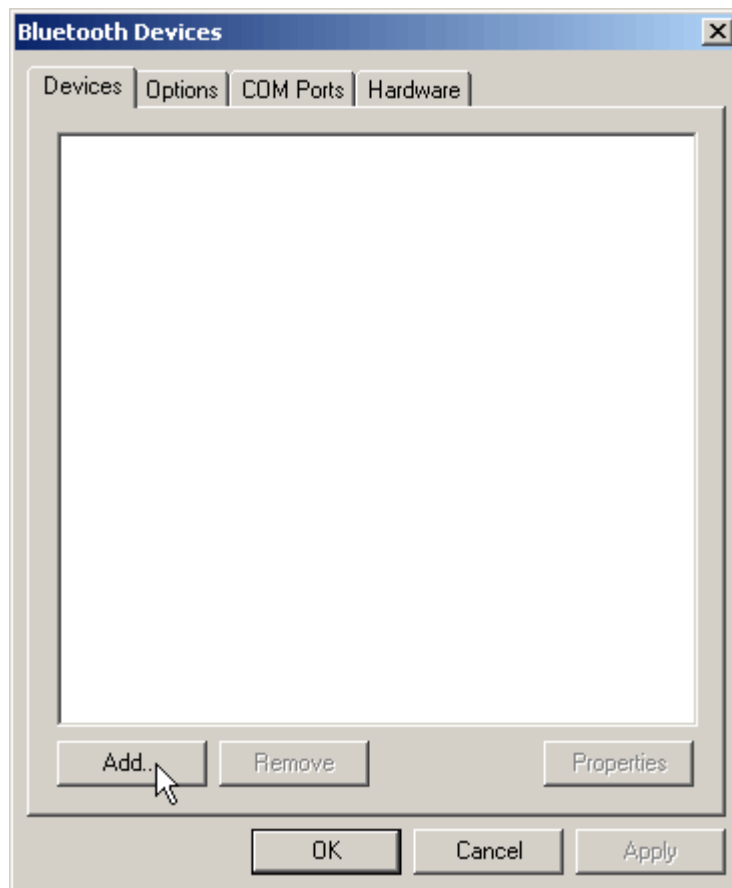

3. The **Add Bluetooth Device Wizard** will open. Make sure that g.MOBIIab+ is switched on and select the checkbox **My device is set up and ready to be found** and click **Next**.



4. The wizard will search for all Bluetooth devices in the range. Select the Bluetooth icon with the serial number of your g.MOBllab+ (e.g. MP-2007_12_01) and click **Next**.

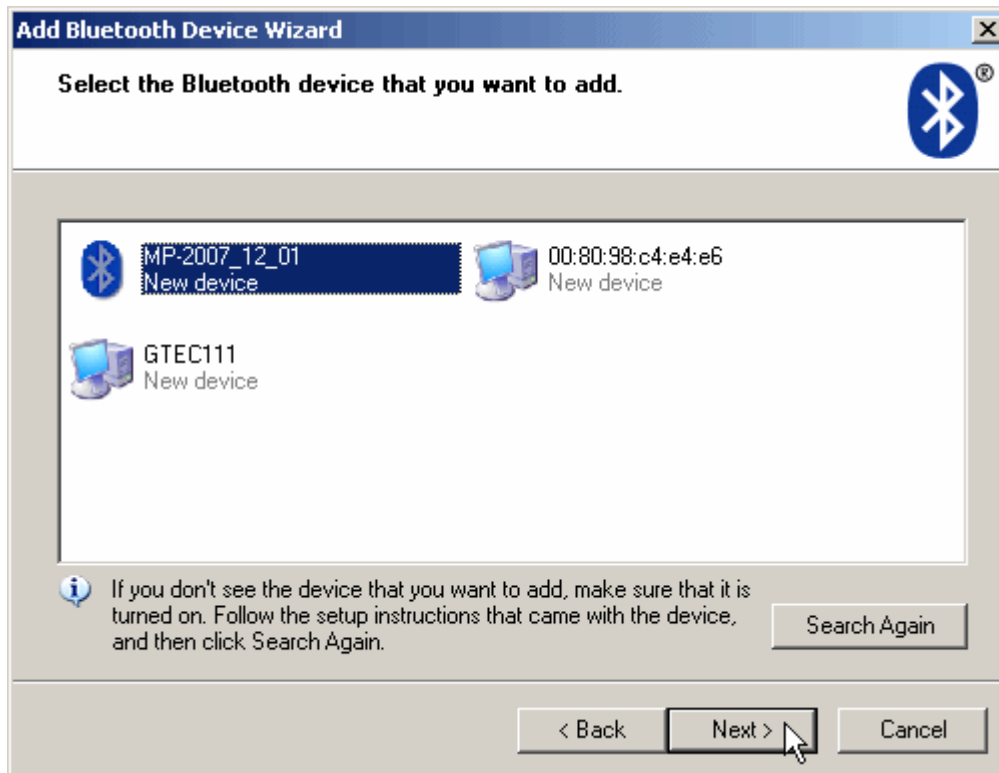

5. A passkey is needed to establish a connection to g.MOBIIab+. Select the radio button **Use the passkey found in the documentation** and enter 1234 as passkey. Click **Next**.

Add Bluetooth Device Wizard

Do you need a passkey to add your device?

To answer this question, refer to the "Bluetooth" section of the documentation that came with your device. If the documentation specifies a passkey, use that one.

☐ Choose a passkey for me

☒ Use the passkey found in the documentation:

☐ Let me choose my own passkey:

☐ Don't use a passkey

You should always use a [passkey](#), unless your device does not support one. We recommend using a passkey that is 8 to 16 digits long. The longer the passkey, the more secure it will be.

< Back **Next >** Cancel

6. The **Outgoing COM port** (e.g. **COM7**) is the port used for every connection to g.MOBllab+. Click **Finish** to exit the wizard.



7. The COM port settings can be reviewed by opening **Bluetooth Devices** in the **Control Panel** in the **Start** menu.

